



The DDS Cervitrac is the most comprehensive Cervical Traction Device on the market today. It's the first device designed to offer the option of pneumatic traction while in the supine position or upright for the freedom of mobile treatment without being restricted to the confines of a door.

The DDS Cervitrac is comfortable, lightweight and easy to use. The soft traction collar can detach from the Cervitrac frame to provide the user with unrestricted ambulatory freedom while at home, at the office or even while traveling.


It works through the use of patented VAP™ Technology (Vertical Air Pressure). When inflated, the advanced longitudinal pneumatic expansion system lifts from the occipital bone and unloads the cervical spine significantly reducing pressure on the intervertebral discs and on pinched nerves. This acts to reduce associated nerve pain helping the cervical spine return to a more natural physiological state and hydration, thus easing pain and assisting in the patient's active rehabilitation.

Features of the Cervitrac™

- **Ambulatory cervical traction treatment**
- **Reduces pressure on nerves**
- **Promotes optimal therapeutic environment allowing discs to rehydrate**
- **Helps reduce or eliminate cervical pain**
- **Reduces muscle tension / Increases range of motion**
- **Provides in-home / at work clinical grade traction treatment**

Treatable Symptoms:

- **Cervical pain syndrome**
- **Pinched nerve**
- **Cervical Radiculopathy**
- **Herniated / bulging disc**
- **Cervicogenic Headache / Chronic Migraine**
- **Degenerative Disc Disease**
- **Whiplash injuries**
- **Joint fixation**

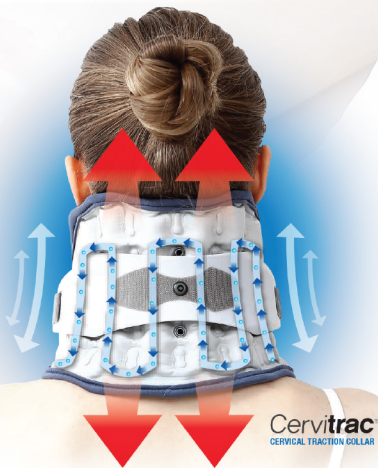

When inflated, our Vertical Air Pressure Technology lifts from the occipital bone and unloads the cervical spine



1 TRACTION COLLAR

The unique and patented expansion air bladder design creates optimal cervical traction.

2 CHIN COLLAR

The contoured frontal collar is soft and ergonomically designed to provide comfort, support and stability.

3 ONE-TOUCH VELCRO TYPE CONNECTORS

Provides simple, easy-to-use fixed-points for quick, simple application and detachment. Stabilizes pressure on the collarbone and chest area.

4 STAND

Provides stability and support for pneumatic traction while in the supine position.

5 DDS HAND AIR PUMP

Standard Manual Hand Air Pump.

**Recommended psi 7-9 lbs
(creating approximately 25 lbs of distractive force)**

SIZE	CIRCUMFERENCE
Small	11.5" - 14.5"
Large	14.6" - 17.5"